


# Complete & Retrofit Companion Care Veterinary X-Ray Systems




Advanced Digital X-ray Solutions for Modern Companion Animal Care.

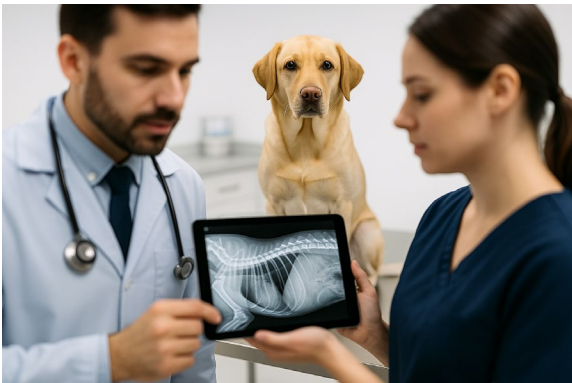

20/20 Imaging delivers digital X-ray solutions designed exclusively for companion animal care. From compact, space-saving systems to high-performance units for busy clinics, our technology helps you capture clear, detailed images for faster, more confident diagnoses. Choose from portable or in-room systems, flexible flat panel detectors, and complete turnkey packages — all powered by intuitive veterinary imaging software built for small animal workflows. With rapid image acquisition, easy patient positioning, and tools to simplify case review, your team can focus more on patients and less on troubleshooting. Backed by reliable performance and expert support, 20/20 Imaging is your partner in providing modern, compassionate care for pets.

## VETSMART™ Table System X-Ray Solution

Veterinary X-ray Table with High-Definition Flat Panel DR — Integrated Camera, Precision Measurement, and Auto APR

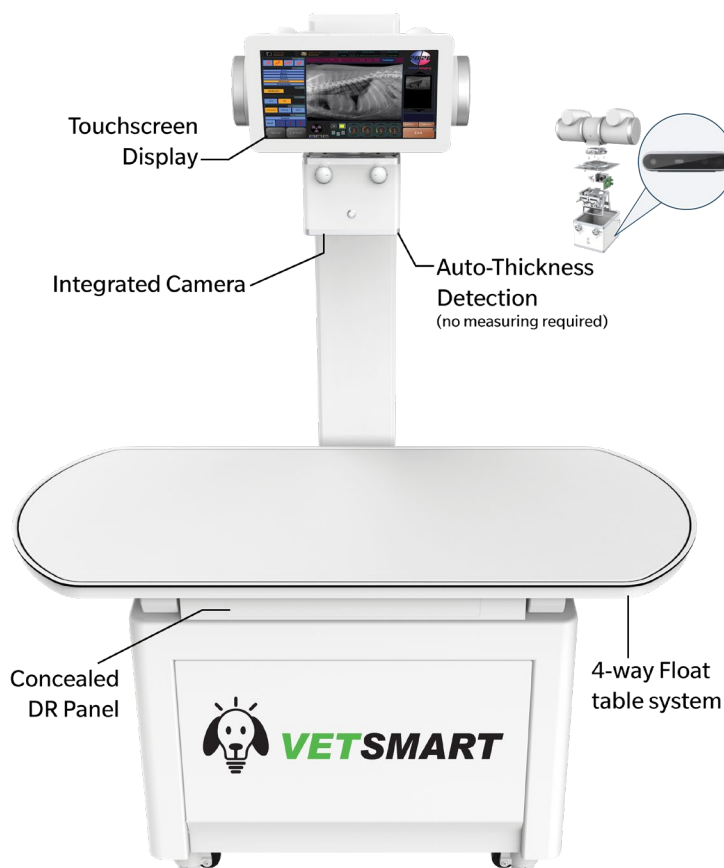


Elevate your veterinary imaging with the VETSMART Digital X-ray System from 20/20 Imaging — a complete, purpose-built solution for companion animal clinics. Designed to combine advanced imaging technology with streamlined usability, VETSMART helps you work faster, smarter, and with greater confidence.

The SMART X-Ray Collimator features patented Auto Thickness Detection (ATD), automatically measuring patient anatomy and optimizing exposure settings for consistent, high-quality images while minimizing radiation dose. Its integrated camera works seamlessly with measurement tools, allowing the tech to view the animal on the table and simply tap where they want to measure, reducing guesswork and saving time.

Available in customizable table configurations — including urine track options and multiple tabletop sizes — VETSMART accommodates a wide range of patients and procedures. Its bright 15" touchscreen makes operation intuitive, while advanced acquisition software delivers sharp, detailed images in seconds. Built-in features like auto-contrast, image auto-shutter/crop, and instant availability enhance diagnostic accuracy and workflow efficiency. Choose the UniversalDR tethered panel or the Momentum DR wireless panel to suit your clinic's needs, both offering exceptional resolution and rapid image capture. Backed by 20/20 Imaging's commitment to quality and support, VETSMART is a smart investment in the efficiency and excellence of your companion animal practice.




### System Dimensions

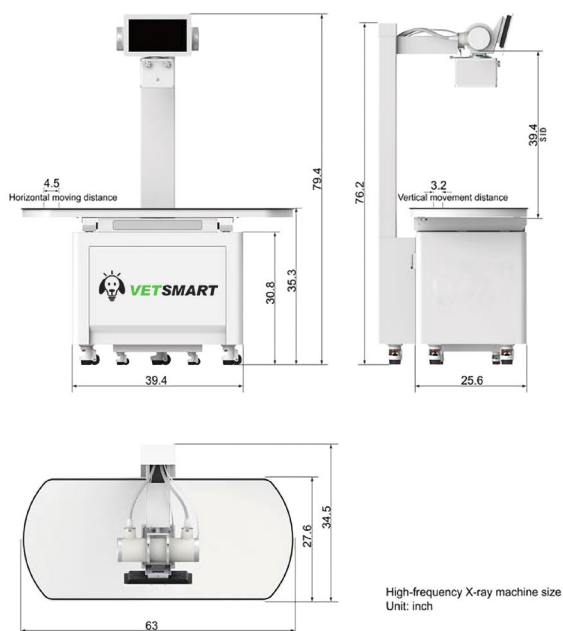


### DR Panel Options

The VETSMART system offers two detector choices: the UniversalDR tethered panel for reliable, fixed connection imaging, and the Momentum DR wireless panel for lightweight portability. Both deliver sharp, high-resolution images in seconds to support fast, confident diagnoses.

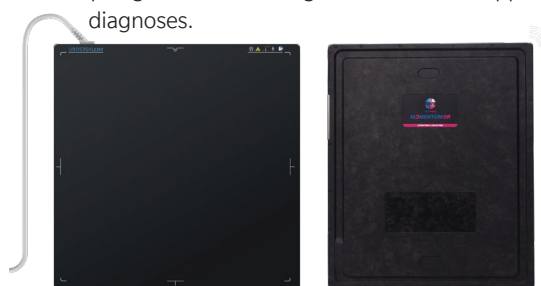


## Digital Retrofit DR Solutions

Digital Imaging Solutions Tailored for Veterinary Care

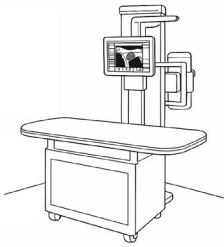

20/20 Imaging offers cost-effective veterinary retrofit solutions that **convert your existing analog X-ray equipment into a fully digital system**. Our retrofit packages integrate seamlessly with most veterinary X-ray units, letting you upgrade to high-resolution digital imaging without replacing your current hardware. Combined with our intuitive software platform for image acquisition and PACS, you'll enjoy faster workflow, enhanced image quality, and veterinary-specific imaging tools — all while maximizing the value of your current investment. It's the smart, budget-friendly way to go digital.

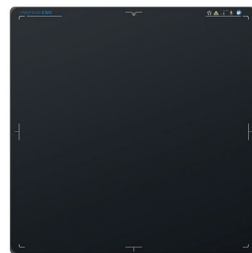

### UniversalDR 17x17 Tethered—Wired (left)

- Always powered through direct connection — continuous operation
- Simple, reliable design for durability and consistent performance
- Scintillator: Csl, 139 micron, 8.8lbs
- 5-Year Warranty

### UniversalDR 17x17—Wireless (Tether-Ready) (center)

- Dual-mode flexibility: wireless detector with optional tether backup
- Balanced mix of wireless convenience and wired reliability
- Scintillator: Csl, 139 micron, 7.5 lbs
- 5-Year Warranty
- Also available in 14x17

### Momentum DR GL-P 14x17—Wireless (right)

- Lightweight Glassless design for durability and reliability
- Built-in capacitor, fully charges in 15 minutes — no batteries to replace
- REALISM™ image processing
- Scintillator: Csl, 200 micron, 3.9 lbs
- 3-Year Warranty

◆ Up to 7 Years Extended Panel Warranties Available\*


\*UniversalDR 17x17 panels include a 5-year warranty, extendable to 7 years. Momentum DR 17x17 panels include a 3-year warranty, extendable to 7 years.

## Your Complete Veterinary Imaging Software Solution

### Opal-RAD Veterinary Imaging Software

Powered by Opal-RAD, our veterinary imaging solution offers a fast, intuitive platform for image acquisition, review, and diagnosis. Designed with the needs of companion animal practices in mind, it delivers advanced viewing tools, seamless digital X-ray integration, and veterinary-specific measurement capabilities.

- Fully customizable user and group settings
- Touchscreen-friendly, user-friendly interface
- Import and export x-ray images, custom query & search
- DICOM send/receive — seamless integration with other imaging systems

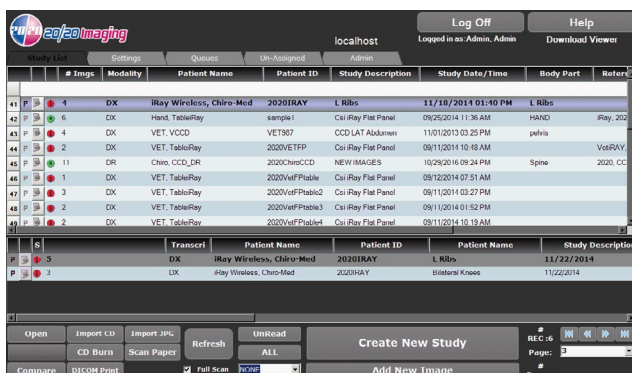

### Opal-RAD Veterinary Image Viewer



Opal-RAD's veterinary image viewer delivers clear, detailed imaging with tools designed for small animal diagnostics. Veterinarians can adjust brightness, contrast, zoom, and orientation, plus apply measurements and annotations directly on-screen. A dedicated **Vertebral Heart Size (VHS) tool** streamlines cardiac evaluations, enabling quick, accurate heart size assessments.



## Acquisition Software for Veterinary Imaging



Our Acquisition Software delivers a seamless, powerful platform for veterinary imaging — whether paired with a complete X-ray system, a digital panel upgrade, or a retrofit solution. Designed to simplify workflow from capture to review, it offers intuitive tools, smart

automation, and full integration with veterinary imaging needs.

- Under 6-second image processing for faster patient turnaround
- Complete DR console for system control and image acquisition
- Automated image processing for fast, diagnostic-quality results
- Customizable technique settings for static or motion studies
- Integrated positioning guides for consistent, high-quality imaging

## Also Available – Portable Vet Solutions



For equine, mixed-animal, and mobile veterinary imaging, the Momentum-EQ (10×12) and Momentum-MA (14×17) systems deliver fast, high-quality digital X-ray in versatile floor-mounted or mobile designs. Lightweight and all-in-one, they're built to handle everything from large equine cases to everyday small animal care.


## mOpal—Zero-Footprint Web Viewer

Access and view studies anytime, anywhere with our mobile-friendly web viewer. Designed for speed and simplicity, the platform runs on any smart device—no downloads or plugins required.

- Works on iPad®, iPhone®, Android™, tvOS & more
- Lightweight, fast, and intuitive interface
- Tap, drag, pinch to navigate
- Patient filtering, search, and basic annotations
- Instant access to viewer and study list



## REALISM™ Image Processing

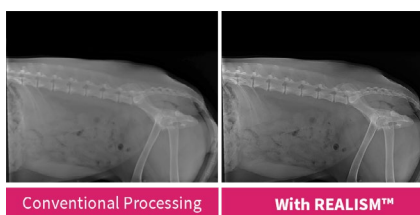

The Momentum DR wireless detector's advanced sensitivity lowers the required X-ray exposure, enhancing safety without compromising image quality. REALISM

technology sharpens high-frequency microstructures—such as bone trabeculae—ensuring unparalleled clarity. Combined with the panel's high resolution, even the finest details are brought to life for more confident diagnoses.

20/20 Imaging is a leading provider of **digital imaging solutions** tailored for veterinary care.

With a full-service approach—from advanced hardware to innovative software—20/20 Imaging streamlines workflows and enhances diagnostics. As part of Konica Minolta, trusted widely by hospitals, we leverage the same proven software, bringing exceptional value to your practice.

Learn more about 20/20 Imaging's **digital x-ray solutions** at [2020imaging.net](https://2020imaging.net).



Learn more about our veterinary solutions at  
[2020imaging.net/vet](https://2020imaging.net/vet)

© 2025 Konica Minolta Healthcare Americas, Inc. M2390 0925 RevA

20/20 Imaging, a division of Konica Minolta Healthcare Americas, Inc.  
829 S Virginia Rd - Suite A  
Crystal Lake, IL 60014  
Tel: +1 (866) 734-6234  
[info@2020imaging.net](mailto:info@2020imaging.net)

[2020imaging.net](https://2020imaging.net)